

Style* Follow

SKINABLE
GALVANIC
Manual

Style Follow Galvanic is a premium skin care device that removes dead skin and wastes without irritation to the skin, helps to absorb deeply the nutrients of cosmetics that are used in basic cosmetics, and makes skin moist and elastic with a micro vibration that vibrates about 30,000 times per second.

STEP1 It penetrates deeply into the skin using shock waves generated when ultrasonic waves are propagated in water, and removes waste materials on the surface of keratin and sebum and remnants after cleansing

Peeling & Sebum & Deep cleansing

1 time of use 3 minutes, recommended once to twice a week

1. If you press the button on the product for 2-3 seconds, the peeling mode is activated
2. Use after cleansing and with moisture remaining on the skin
3. Lightly touch the tip of the blade to the skin, then use it as if you are pulling moisture away from outside the skin

STEP2 Galvanic ions and vibration function calms sensitive skin and absorbs active ingredients more deeply, further enhancing effectiveness

Nutrition Care Soothing & moisturization

1 time of use 3 minutes, available daily

1. Pressing the button in peeling mode will activate the nutrition mode after about 1 seconds
2. Use with cosmetics (ampoule or essence) on skin
3. Apply a large portion of the blade lightly to the skin, then slowly massage out of the skin texture

STEP3 Galvanic ion and vibration function absorbs active ingredients, stimulates skin tissue and pores, helps to increase elasticity

Tightening Care Wrinkle care & Rejuvenate


1 time of use 3 minutes, available daily

1. If you press the button in the nutrition mode the tightening mode is activated after about 1 second
2. It is used for areas where elasticity care is needed after the nutrition mode
3. Apply a large portion of the blade lightly to the skin then slowly massage out of the skin texture.

STEP1



STEP2

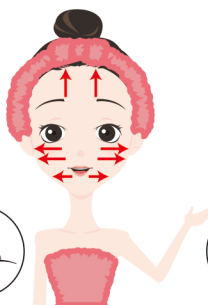

STEP3



How to Operate

Power ON/OFF

- * The Green LED lights up when the button is pressed. .
- * If you press the button for about 2-3 seconds, it will activate.
- * Press the button for 2 ~ 3 seconds to turn off the power.

Switching Mode

- * The mode changes after about 1 seconds when the button on the product is pressed after operation starts

Auto Shutdown

- * Each mode will automatically turn off after 3minutes


How to Charge

- * It is recommended that you charge the product fully after you have received the initial product
- * When the charging cable is connected to the bottom of the product, the LED lights up when charging
- * The product does not work during charging
- * If there is no battery during use, it will automatically shut down

Manage/Storage

- * After use, wipe gently with a dry towel.
- * Never wash with water. It may cause malfunction.
- * Do not expose it to direct sunlight or store it in a humid place.
- * Do not store together with sharp objects that may cause damage.

1. Blade
2. Hand Electrode
3. LED Mode
4. Power On/Off & Mode Switching button
2times - Nutrition mode / 3times - Tightening mode
When using, touch the hand electrode

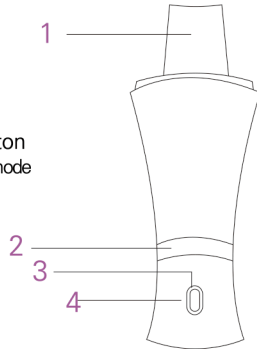

[Precautions]

1. Be sure to read and understand the precise usage before use.
2. Do not use if there is any injury or wound on the skin.
3. If any abnormality occurs, discontinue use and contact the customer center or the seller.
4. To reduce the risk of shock, fire or other injury, do not immerse the appliance in water.
5. Do not place or store the product near bathtubs, showers, washbasins, or other bathrooms.
6. Do not operate the product if it is damaged.
7. Do not immerse in cosmetics because it is not a waterproof product.
8. Please keep out of reach of children.

[Warning]

1. Please stop using when overheating during use.
2. Do not use if the power cord is damaged.
3. Do not use it continuously for a long time.
4. Avoid use when you have received heart surgery, heart disease, fever, infectious disease, cancer, cuts, acute illness, hypotension, cardiac pacemaker, artificial heart or lungs, skin disease, burns, or metal implants or implants..

[Consumer Error]

1. Failure due to user error or carelessness (falling, flooding, shock, breakage, unreasonable operation, etc.).
2. In the event of malfunction due to consumer intention or negligence.
3. In the event of malfunction or damage to the product due to use of non-genuine parts or accessories.
4. Discolored by leaving the cosmetics left on the product.

Style Follow Limited Warranty

Warranty Card			
Product	Galvanic Massager	Model Name	Skinable Galvanic
Purchase Date		Warranty Period	One (1) Year
Customer Name		Contact Number	
Customer Add			

What Does This Limited Warranty Cover?

This limited warranty covers any defects in material or workmanship in the new Skinable Product accompanied by this limited warranty statement. Only consumers purchasing this product from an authorized retailer or reseller may obtain coverage under this limited warranty. The warranty period is one(1) year from purchase date. (eg. by presenting proof of purchase).

What Does This Limited Warranty Not Cover?

This limited warranty does not cover any problem that is caused by (a) commercial use, accident, abuse, neglect, shock, electrostatic discharge, degaussing, heat or humidity beyond product specifications, improper installation, operation, maintenance or modification; or (b) any misuse contrary to the instructions in the user manual; or (c) malfunctions caused by other equipment. This limited warranty is void if a product is returned with removed, damaged or tampered labels or any alterations (including the unauthorized removal of any component or external cover). This limited warranty does not cover, also, consequential damages; incidental damages; and costs related to data recovery, removal, and installation are not recoverable under this warranty. This limited warranty applies only to hardware products that are licensed (and warranted where applicable) pursuant to separate written agreement

IMPORTANT RESTRICTION

Skinable may restrict warranty service to the country where Skinable or its Authorized Distributors originally sold the Skinable Product.

YOUR RESPONSIBILITIES

Before receiving warranty service, Skinable or its agents may require that you furnish proof of purchase details, respond to questions designed to assist with diagnosing potential issues and follow Skinable procedures for obtaining warranty service.

Following warranty service your Skinable Product or a replacement product will be returned to you as your Skinable Product was configured when originally purchased.

Important: Do not open the Skinable Product. Opening the Skinable Product may cause damage that is not covered by this warranty. Only Skinable should perform service on this Skinable Product.

GENERAL

No Skinable reseller, agent, or employee is authorized to make any modification, extension, or addition to this warranty. If any term is held to be illegal or unenforceable, the legality or enforceability of the remaining terms shall not be affected or impaired. This warranty is governed by and construed under the laws of the country in which the Skinable Product purchase took place. Skinable is identified at the end of this document according to the country or region in which the Skinable Product purchase took place. Skinable or its successor in title is the warrantor under this warranty.